


How to improve the reproducibility of your flow chemistry

To enhance the reproducibility of flow chemistry scale-up - precise control of mixing and temperature in highly exothermic or rapid reactions is essential.

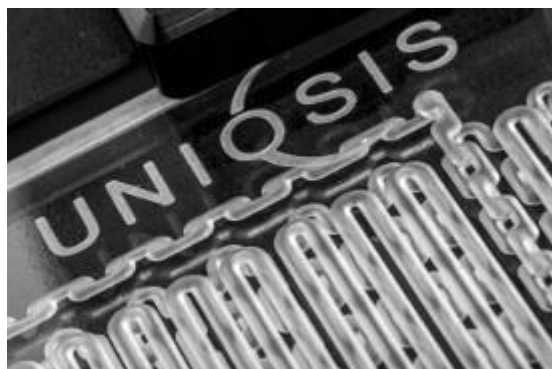

Glass Static Mixer (GSM) chip reactor blocks

The expanding range of **Glass Static Mixer (GSM) chip reactor blocks** from **Uniqsis** produce an efficient, turbulent mixed reagent stream for flow chemistry reactions that is not diffusion dependent.

Available in sizes

from 270 μ l to 20 ml, the Uniqsis range of high-quality GSM chips can be used to perform reactions from -80°C to over 150°C. GSM chips up to 2 ml can be operated up to 40-bar and are available in 2-channel or 3-channel inlet configurations. Larger GSM chips of 10ml and 20 ml have a premixing channel followed by a residence domain and can operate up to 10-bar.

Glass Static Mixer chips

can be used as mixing modules prior to a coil reactor residence time unit, or as reactor blocks for rapid exothermic reactions.

GSM chips

from Uniqsis are precision machined from inert borosilicate glass to withstand a wide temperature range and can be conveniently attached to a FlowSyn column heater module, Cold Coil or Polar Bear Plus cryogenic reactor module.

For further information

please visit https://www.uniqsis.com/paProductsDetail.aspx?ID=ACC_CHIP or contact Uniqsis on +44-1223-942004 / info@uniqsis.com



Uniqsis Ltd.

Since 2007, has specialised in the design and supply of mesoscale continuous flow chemistry systems for a wide range of applications in chemical and pharmaceutical research. The company's aim is to make flow chemistry easily accessible to both novices and experienced users.

Global HQ

Uniqsis Ltd

Foxes Bridge Farm Unit 1,
Lime Tree Barn
Foxes Bridge Farm
Royston Lane
Comberton
Cambridgeshire CB23 7EE
UK

tel: +44-1223 942004

email: info@uniqsis.com

web: www.uniqsis.com